NOCOSPRAY

Kube

- FR** Manuel d'Utilisation
- EN** User Manual
- ES** Manual de instrucciones
- DE** Gebrauchsanweisung
- IT** Manuale d'uso
- NL** Handleiding

FR

SOMMAIRE

FR

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4.2	Device positioning.....	41
4.3	Filling the bottle	41
4.4	Starting the diffusion process	42
4.5	LED indicator light.....	42
5.	TRACKING SYSTEM	43
5.1	Report export.....	43
5.2	Report template.....	44
6.	POWER MANAGEMENT	44
7.	APPENDICES	45
7.1	Maintenance	45
7.2	Technical specifications	46
7.3	Error messages.....	46
8.	RECYCLING	47
9.	NOTES	47

1. INTRODUCTION

EN

1.1 General information

NOCOSPRAY KUBE is a portable device used for the diffusion of disinfectant products.

This device is for professional use only. It can be used in hospitals, the food industry, the pharmaceutical industry, or in any environment where disinfection is required.

The NOCOSPRAY KUBE operates on battery power only.

This user manual describes the operation of the NOCOSPRAY KUBE diffuser and includes essential information, guaranteeing safe use of the device.

This documentation regards NOCOSPRAY KUBE diffusers in their original state, without any modifications.

- Carefully read this manual before using the NOCOSPRAY KUBE diffuser. It has been created to allow you to understand and efficiently use all of the device's features.
- The illustrations in this manual are there to help you understand the various functions of the NOCOSPRAY KUBE. They are not binding.

1.2 Classification / Regulations

The NOCOSPRAY KUBE diffuser complies with the following standards and directives:

Electromagnetic Compatibility (EMC): Directive 89/336/EEC Standard EN61000-4-3 (2006 + A2:2011)	The product does not emit any electromagnetic disturbances which may affect other devices. It presents appropriate immunity levels in the presence of electromagnetic disturbances.
Electrical safety: Low Voltage Directive 2014/35/EU Standard EN 61010-1 (2011 + A1:2019)	When used correctly and in compliance with its intended use, the product does not compromise user health and safety.

1.3 Guarantee

The guarantee certificate must be presented to the retailer or to OXY'PHARM - Z.A. des Grands Godets – 829, rue Marcel Paul – 94508 Champigny-sur-Marne, France, to prove the validity of the guarantee.

Guarantee clauses

The following provisions are not exclusive of the entitlement to the legal guarantee provided for under Articles 1641 et seq. of the French Civil Code pertaining to hidden defects.

Purpose of the guarantee

The guarantee consists of the free provision of parts identified as being defective by our technical departments. Repairs made following improper or abnormal use, negligence or overload of the device, or those required due to variations in the electrical network, power surges or defective installations, etc., are not covered under the guarantee.

In addition, if abnormal corrosion of soft parts, hoses or seals, etc., is identified, signifying the use of a liquid other than those produced by OXY'PHARM, the guarantee will be declared null and void.

Transport

In the event of damage during transport, the recipient must report all reservations regarding the carrier before taking delivery of the device.

The device must be stored in a dry place and transported in its original packaging.

Guarantee period

The NOCOSPRAY KUBE is guaranteed for a 2-year period from the invoice date.

2. SAFETY GUIDELINES

The manufacturer may not be held liable for any material or bodily damage caused by improper use of the device in contravention of the precautions recommended in the non-exhaustive list below:

	<p><i>The device has a vibrating nozzle. Its membrane must not be touched.</i></p> <p><i>Do not insert any objects or liquids into the nozzle, either while the device is in operation or switched off.</i></p>
	<p><i>To avoid any risk of electrical shock, while it is charging, the device must be connected to an electrical supply network which is fitted with a protective earth connection.</i></p> <p><i>Make sure that the electrical supply complies with the applicable standards.</i></p>
	<p><i>The device must never be used when it is connected to the mains power.</i></p>
	<p><i>This device is intended for professional use only.</i></p> <p><i>Any modifications made to the EM device are prohibited.</i></p> <p><i>Portable and mobile RF communication devices may affect the operation of EM devices.</i></p> <p><i>Using accessories, transducers and cables other than those specified, with the exception of equipment sold by the manufacturer as replacement parts for internal components, may lead to an increase in emissions or a reduction in the device's immunity levels.</i></p> <p><i>The device must not be used next to or stacked with other devices</i></p> <p><i>If this condition cannot be respected, test that the device operates properly in this configuration.</i></p> <p><i>Make sure that the electrical supply complies with the usual standards.</i></p>
	<p><i>Only use the product provided by Oxy'Pharm.</i></p>
	<p><i>Remove the bottle from the diffuser when transporting the device.</i></p> <p><i>Handle the container with care.</i></p>
	<p><i>Users are not authorised to carry out repairs on the device.</i></p> <p><i>If the device stops working, call your retailer/approved technician.</i></p>
	<p><i>Any electrical work on the device must be carried out by an approved technician. OXY'PHARM provides retailers with all of the technical information necessary for troubleshooting.</i></p>

ENTER THE CONTACT DETAILS OF YOUR RETAILER HERE

Name:

Address


Tel.

The device diffuses a chemical substance. Given this, the following recommendations must be followed:

	<i>All accompanying documents must be read. The decontamination process must be monitored.</i>
	<i>The device must only be used indoors.</i>
	<i>Do not inhale the product.</i>
	<i>Do not smoke close to the device.</i>
	<i>Do not direct the device towards an incandescent zone.</i>
	<i>Do not use in the presence of a flammable anaesthetic mixture.</i>
	<p><i>The decontamination process must be monitored. Leave the room and close the door while the machine is in operation, and during the contact time, and do not re-enter the room until after these periods. There must be no humans present during the treatment process.</i></p> <p><i>Do not obstruct the ventilator grill.</i></p> <p><i>Respect the contact times recommended in the protocols for use.</i></p> <p><i>Respect a minimum 30-minute period after completion of treatment before re-entering the treated room.</i></p> <p><i>In case of emergency, if someone has to enter the treated room, they must wear Personal Protection Equipment: In order to ensure the safety of personnel who are not protected by PPE, a H2O2 detector may be used to measure any residual gas. The room may only be entered when the H2O2 concentration (ppm) is below the required levels stipulated by local laws regarding the authorised exposure times in professional settings.</i></p>

3. DESCRIPTION

3.1 Description of the NOCOSPRAY KUBE

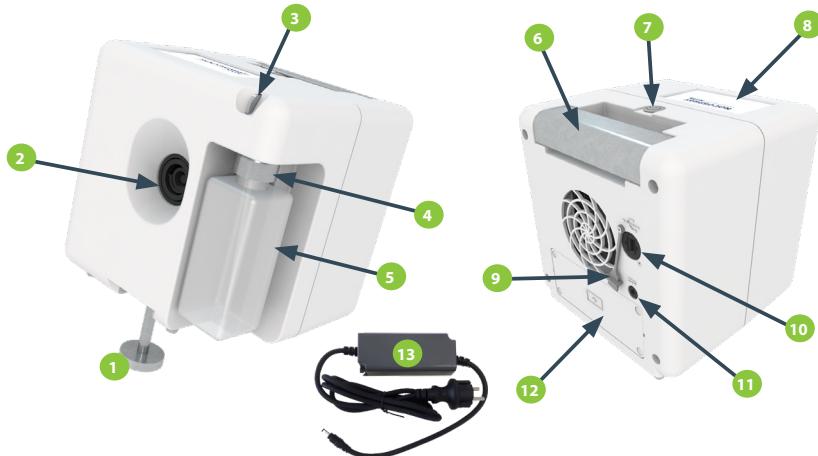1 Height-adjustable leg

2 Diffusion membrane

3 Bottle unlock button

4 Cap

5 Bottle

6 Handle

7 General illuminated switch

8 Control touch screen

9 USB port cover

10 USB port


11 Power cable plug

12 Battery compartment

13 Charger

3.2 Dimensions

EN

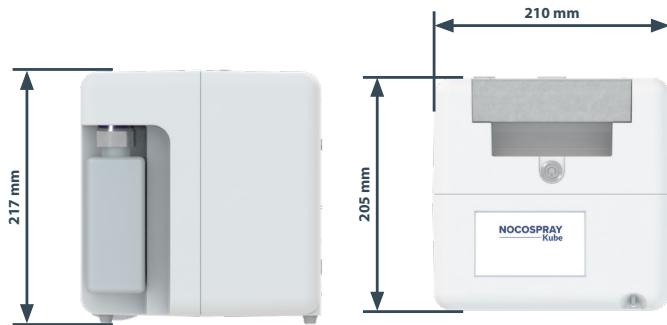4. USE

4.1 Device control panels

4.1.1 Home screen

The home screen lets you programme the device and view the settings entered. The time, battery level, device serial number, and software version are all displayed at the top of the screen.

Illustration	Action
	"Operator" menu The selected user is displayed next to the icon.
	"Location" menu. The selected location of protocol execution is displayed next to the icon.

Illustration	Action
	"Configuration" menu for the diffusion.
	Diffusion settings display: <ul style="list-style-type: none">• Volume to disinfect.• Diffused dosage.• Delay time.• Contact time.
	Programme start.
	"Settings" menu.

4.1.2 "Operator" menu

The "Operator" screen lets you enter the name of the operator and the company name.

The operator information is not mandatory. This information is saved to the device's memory.
In locked administrator mode, you cannot add or remove an entry.

4.1.2.1 "Operator" menu, general screen

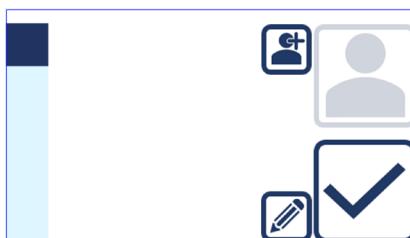
Illustration	Action
	Add an operator
	Edit the operator selected from the list.
	Selection confirmation.

4.1.2.2 "Operator" menu, edit screen

EN

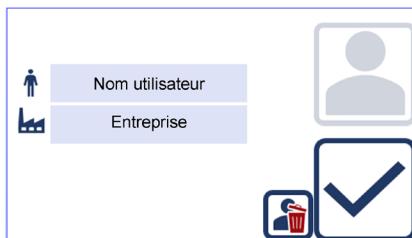
Illustration	Action
	Remove an operator
	Selection confirmation.

4.1.3 "Location" menu

The "Location" screen lets you enter the name of the lab and the unit name.

The location information is not mandatory. This information is saved to the device's memory.
In locked administrator mode, you cannot add or remove an entry.

4.1.3.1 "Location" menu, general screen

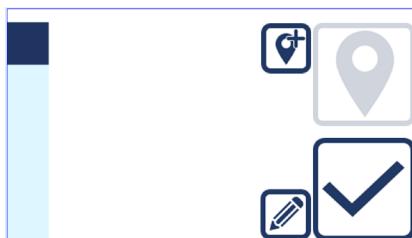
Illustration	Action
	Add a location
	Edit the location selected from the list.
	Selection confirmation.

4.1.3.2 "Location" menu, edit screen

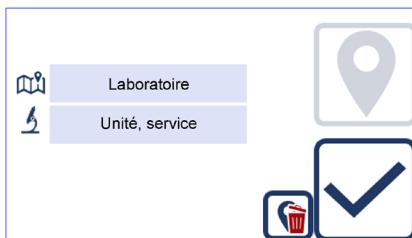
Illustration	Action
	Remove an operator
	Selection confirmation.

4.1.4 "Configuration" menu

The "Configuration" screen lets you enter the name of the location, serial number, volume, dosage, cycle number, and contact time.

The Configuration information is mandatory. This information is saved to the device's memory.
In locked administrator mode, you cannot add or remove an entry.

4.1.4.1 "Configuration" menu, general screen

EN

Illustration	Action
	Add a configuration.
	Edit the configuration selected from the list.
	Selection confirmation.

4.1.4.2 "Configuration" menu, edit screen, page 1

Illustration	Action
	Name of the zone or equipment to be disinfected.
	Serial/ID number of the zone or equipment.
	Volume to disinfect (from 0.1 to 15 m³).
	Next page.

4.1.4.3 "Configuration" menu, edit screen, page 2

EN

Illustration	Action
	Previous page.
0 ml/m ³ x 0 cycles	Product dosage and number of cycles.
	Interphase time between each cycle.
	Delay time before diffusion start.
	Contact time.
	Remove the configuration.
	Selection confirmation.

4.1.5 “Settings” menu

The “Settings” screen lets you enter/edit the device’s information.

4.1.5.1 “Settings” menu, general screen

Illustration	Action
*****	Creation/change of administrator password.
<input type="checkbox"/>	Unlock administrator mode.
08:19:30	Device time.
21/06/2022	Configure/edit the device time and date.
094538192	Device serial number.
521 min	Remaining duration of device battery in use
Français	Device language.
	Remove reports from the device. If an administrator password is entered, this function is password-protected.
	Bottle replacement procedure (enter the product information: Batch number, expiry date, etc.).
<input type="checkbox"/>	Enabling the remote control (option disabled).
<input type="checkbox"/>	Enabling temperature and humidity recording (allows the operator to indicate these values just before starting the cycle via external equipment, which are under their responsibility).

4.1.5.2 "Settings" menu, time and date screen

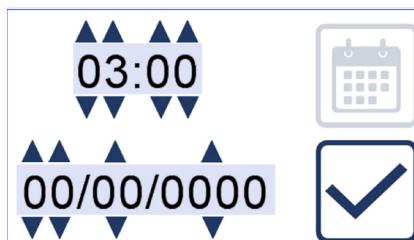
Illustration	Action
	Use the arrows to edit the value.
	Selection confirmation.

4.1.5.3 "Settings" menu, "Product information" screen

Illustration	Action
	Enter the product name.
	Enter the product batch number.
	Enter the product expiry date.
	Cancel the information.
	Confirm the information.

If the bottle isn't replaced and is just manually refilled instead, it must be fully refilled.
The NOCOSPRAY KUBE system registers a full refill as being equivalent to 225 ml (i.e., 250 ml -10% safety margin).

4.2 Device positioning

Treating a room:

The device must be positioned at mid-height in the room. The diffusion must be directed towards the centre of the room.

Treating equipment:

The device must be positioned at the base of the equipment. The diffusion must be directed upwards, at a diagonal angle.

The NOCOSPRAY KUBE must be positioned on a flat and resistance surface.

Only carry the device by its handle.

4.3 Filling the bottle

To remove and insert the bottle, follow the steps shown below.

When changing the bottle, remember to enter the product name, batch number and expiry date by pressing the icon in the "Settings" menu

Always check that the bottle contains sufficient product for the treatment to be carried out.

Always use a full bottle to replace it, if needed.

Wear appropriate gloves when replacing the bottle.

In the event of a product spill, wipe up the spilled product immediately and place the paper towels used in a container/basin/bucket containing water in order to neutralise the chemical reaction.

Steps	Action	Illustration
1	Release the bottle cap by pressing the button and pulling downwards.	
2	Unscrew the cap from the bottle in the device.	
3	Date the new bottle of product, then screw the cap back on.	
4	Put the bottle back into the device. You should hear a "click".	

4.4 Starting the diffusion process

The device must never be used when it is connected to the mains power.

When the configuration for the diffusion process is complete, the information will appear on the home screen.

Steps	Action	Illustration
1	Confirm the process.	
2	The safety screen will be displayed.	
3	Start the process. The delay time is requested again before the diffusion starts.	
4	The diffusion cycle, interphase, and end times are all displayed.	
5	At the end of the process, the treatment end screen will be displayed.	

If needed, the device may be stopped during the treatment process by pressing

4.5 LED indicator light

Steps	Process
1	Delay time: Red flash (50 ms red, 500 ms off)
2	Diffusion in progress: Blinking red (500 ms red, 500 ms off)
3	Contact time: Blinking red (500 ms red, 500 ms off)

Steps	Process
4 (a)	End of process, if compliant: Solid green light until you go back to the home screen
4 (b)	End of process, if non-compliant: Solid red light until you go back to the home screen

5. TRACKING SYSTEM

5.1 Report export

Steps	Action	Illustration
1	Insert a USB stick into the device's USB port, after removing its cover.	
2	The report export window opens automatically.	
3	Confirm the extraction.	
4	Remove the USB stick once the export is completed, and remember to put the cover back over the port.	

Elements are classed on the USB stick as follows:

- "KUBE_DATE_OF_DIFFUSION_NUMBER" folder
- "KR_DATE OF DIFFUSION_TIME OF DIFFUSION" file.



5.2 Report template

EN

08/01/2001 22:43


012345678

NOCOSPRAY Kube

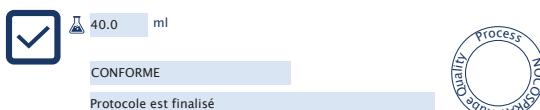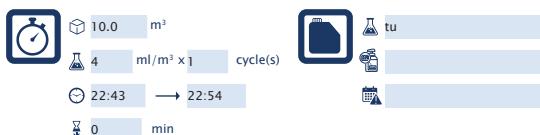6. POWER MANAGEMENT

When the battery level is low, the following screen will suggest that you plug the device into the mains power supply:

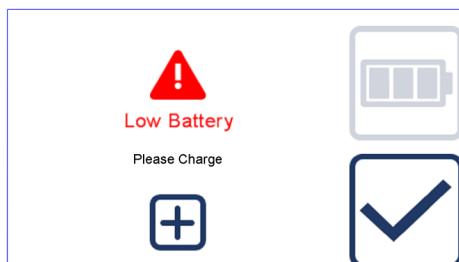If the device stops due to an empty battery:

- When the device is connected to the mains power, it will switch on and display the last screen displayed prior to its shutdown, and the battery will be charged.

End-of-process protocol:

- LED and screen stay illuminated for 2 hours.
- After 2 hours, the screen switches off but the LED light will remain lit until the end-of-process has been validated.

Energy conservation after the end-of-process message display:

- The screen switches off every 5 minutes if the operator doesn't validate the end of the process.
- When it's switched back on, the LED and screen display the last status pre-shutdown.

*Only use the charger provided with the equipment.
In the event of a battery defect, please contact the Product Support Service and do not replace it yourself.*

The device must never be used when it is connected to the mains power.

7. APPENDICES

7.1 Maintenance

*Before any work is carried out, switch off the device and unplug the mains power cord.
Do not use abrasive products, bleach, acetone or solvents for maintenance or disinfection operations carried out on the device. These aggressive products risk damaging the device's surface materials.*

To clean the external surfaces, use a slightly damp cloth or a wipe.

Do not touch the nozzle when cleaning the exterior of the device.

7.2 Technical specifications

Transport and storage conditions	
Temperature	0 °C to 45 °C
Humidity	10% to 95% RH without condensation
Atmospheric pressure	700 to 1,060 hPa
Usage conditions	
Temperature	0 °C to 40 °C
Humidity	10% to 90% RH without condensation
Atmospheric pressure	800 to 1,060 hPa
General specifications	
Device weight	3.1 kg
Package weight	5.5 kg
Package dimensions	270x270x315 mm
Compatible packaging	25 cl
Minimum volume for treatment	0.1 m ³
Maximum volume for treatment	15 m ³
Average liquid flow rate	5 ml/min. ± 5%
Electrical specifications	
Power supply voltage	230 V mono
Power supply frequency	50 to 60 Hz
Power	60 W
Nominal intensity	3 A
Battery specifications	
Battery type	NimH
Voltage	14.4 V
Electrical charge	4 Ah
Number of recharge cycles	500

In order to preserve the battery's life span, a full recharge must be carried out at least once every 3 months.

7.3 Error messages

Message	Meaning
Top up with the product	The bottle is empty or not filled enough. Only replace the bottle with a full bottle.
Low battery	The battery charge is at an insufficient level.

8. RECYCLING

EN

IMPORTANT INFORMATION FOR ENVIRONMENTALLY-FRIENDLY DISPOSAL

This product is subject to separate collection for electrical and electronic equipment, in compliance with Directive 2002/96/EC (applicable in European Union member states and other European countries with separate collection systems).

This symbol indicates that, at the end of its life span, the device must be disposed of separately from household waste. The device must be taken to a selective collection site for electrical and electronic devices, or taken back to the retailer when purchasing a new device.

At the end of its life span, the user must take the device to the competent collection structures or to a duly authorised company, otherwise the user incurs the risk of sanctions under waste disposal legislations.

The selective collection prior to the recycling, processing and environmentally-friendly disposal of the used device helps to limit environmental disturbances and health impacts, and also promotes the recycling of the materials composing the device.

Waste consumables must be handled in compliance with the local waste management regulations.

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