

# NOCOSPRAY

## LADYBIRD

The light and portable aerial surface bio-disinfection device.



■ ■ Manufactured in France by Oxy'Pharm according to ISO 13485.

**Nocospray ladybird** is a light and portable aerial surface bio-disinfection device. It ensures aerial surface disinfection or disinsectisation when used along with our Nocolyse or Oxyphy ranges.

A dry fog allows a uniform treatment of every cm<sup>2</sup> of surfaces, as a preventive or curative measure.

Treatment of volumes from 10 to 1,000 m<sup>3</sup>.



### EASY TO USE

Easy adjustment for the volume of the room to be treated thanks to its touchpad.



### ECONOMICAL

Very small amount of product required.



### GUARANTEED TIME SAVING

No preparation required before treatment (except for cleaning) or after treatment (no need to wipe or ventilate).

## COMPATIBILITY

- Disinfectant bottles (1 L): Nocolyse, Nocolyse+, Nocolyse Food, Nocolyse One Shot.
- Insecticide bottles (1 L): Oxyphy, Oxyphy+.

## PRECAUTIONS FOR USE

Read all product information before use. Instructions for use are available on request by email: [info@oxypharm.net](mailto:info@oxypharm.net)


## TECHNICAL CHARACTERISTICS

Dimensions (L x W x H cm)	49.4 x 30 x 34.2 cm
Empty weight	5.9 kg
Packaging	
- Dimensions (L x W x H cm)	60 x 42 x 42 cm
- Weight	9 kg
Maximal treatment volume	1,000 m <sup>3</sup>
Turbine speed	22,000 rpm
Output speed of the mist	80 m/s
Average liquid flow rate	1,000 ml/h
Storage conditions	
- Temperature	0°C to 70°C
- Hygrometry	10% to 95%
Power	1000 W
Average particle size	<5 µ

## REFERENCES

DEVICE	REFERENCE
Nocospray Ladybird 230 V	2000.001
Nocospray Ladybird 115 V	2000.011

## THE BIO-DISINFECTION PROCESS

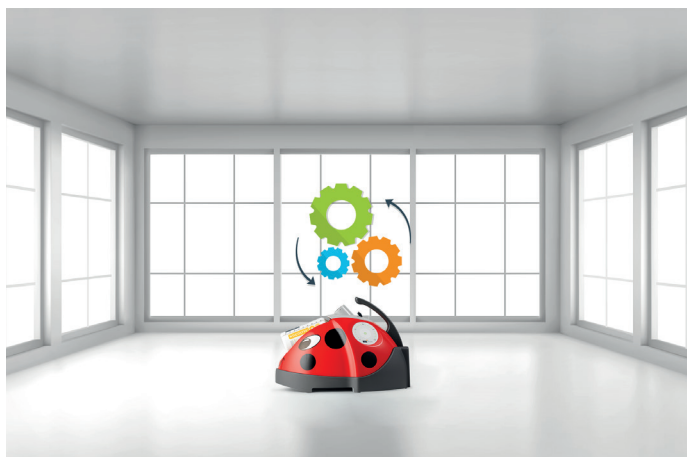


### 1. INSTALLATION

Placement of the Nocospray Ladybird system in the room to be treated.



### 2. SPRAYING

Spraying of the product in the form of a dry mist in the room to be treated, evenly over all accessible surfaces.

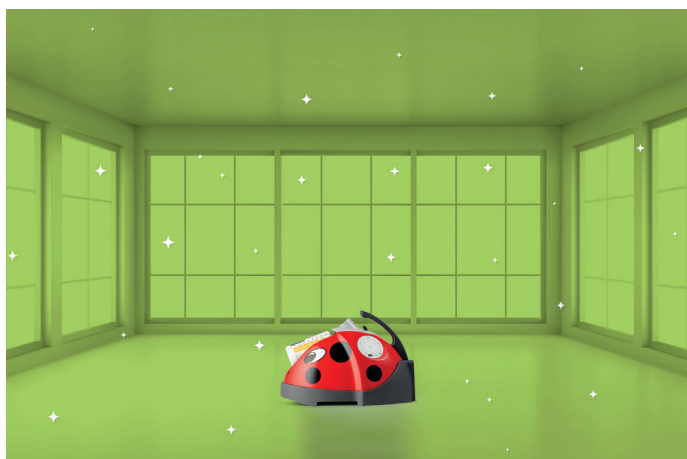

### 3. CONTACT TIME

The dry mist remains suspended in the air and eliminates all micro-organisms (viruses, bacteria, spores, yeasts, etc.) present on the accessible surfaces of the room.



### 4. BIODEGRADATION

The dry mist naturally degrades without leaving any residue and without corroding surfaces (including electronic and mechanical ones). No need to wipe or ventilate.

## ACCESSORIES



#### The programmer

Allows for delayed activation of the device.



#### The wall crossing system

Allows you to treat a room without entering it.




#### Nocotest

Checks the proper distribution of the products on the surfaces.



#### The remote control

For triggering the broadcast from outside the room.

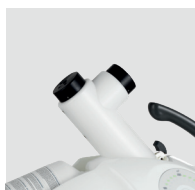

#### The double nozzle

For two-way spraying suitable in long, narrow spaces.



#### Nocobox

Allows control of the microbiological efficiency obtained.